



SHOULDER IMMOBILIZER, PREMIUM

INSTRUCTIONS FOR USE

Single Patient Use Only

BEFORE USING THIS DEVICE, PLEASE READ THE FOLLOWING INSTRUCTIONS COMPLETELY AND CAREFULLY. CORRECT APPLICATION IS VITAL TO THE PROPER FUNCTIONING OF THE DEVICE.

INTENDED USE:

- Shoulder Injuries & Instabilities
- Dislocations

WARNINGS & PRECAUTIONS:

If you experience any pain, swelling, sensation changes, or any unusual reactions while using this product, consult your medical professional immediately.

CLEANING INSTRUCTIONS:

Hand wash with mild soap in lukewarm water, rinse thoroughly. AIR DRY.

Note: If not rinsed thoroughly, residual soap may cause irritation and deteriorate material.


WARRANTY:

Manufacturer guarantees all products to be free from defects in materials and workmanship for a period of 6 months from date of purchase.

PART NUMBERS: 200080-XX

-XX	Size	Length*
-03	SM	Up to 11"
-05	MD	11" – 13"
-07	LG	13" – 15"
-08	XL	15" +

*Measure length from elbow crease to base of index finger

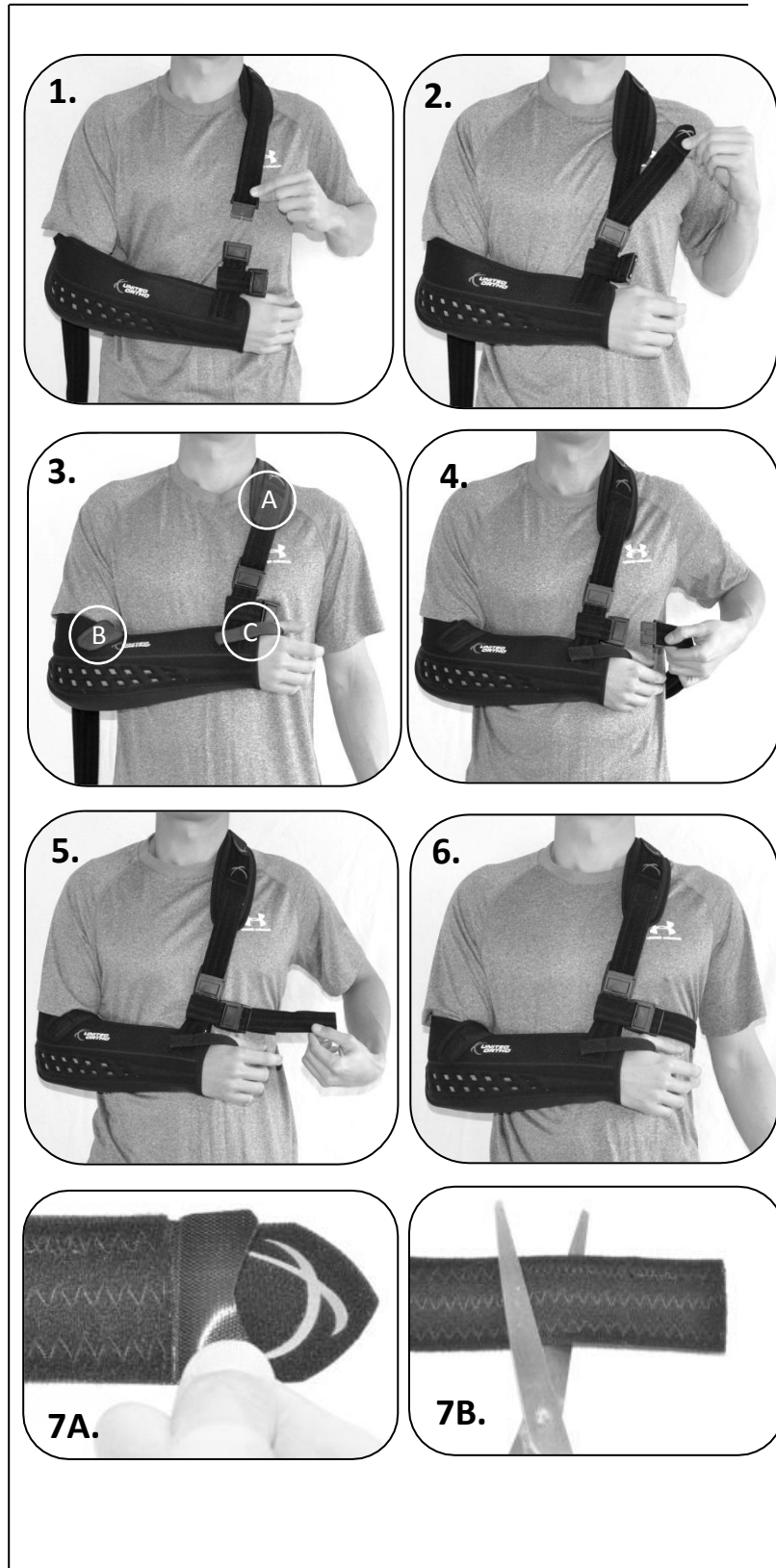

APPLICATION INSTRUCTIONS

1. Detach quick release buckle of shoulder strap. While supporting the injured arm, slide the forearm into the sling. Wrap the adjustable shoulder strap around the neck and shoulders and attach the quick release buckle.
2. The shoulder strap has two adjustable ends. Adjust the shoulder strap to the proper length.
- 3A. The shoulder pad may be adjusted by removing it from the shoulder strap and repositioning it over the neck and shoulder.
- 3B. Attach the closure strap across the top of the sling near the elbow so that it holds the arm securely in the sling.
- 3C. Attach the thumb strap to the open end of the sling between the thumb and fingers.
4. Wrap the waist strap around the back and attach the quick release buckle.
5. The waist strap has two adjustable ends. Adjust the waist strap to the proper length.
6. This is an image of a properly fit brace.
- 7A. Straps may be shortened if too long. Remove the "alligator tab" from the end of the shoulder and/or waist strap.
- 7B. Trim to desired length using scissors. Apply as described above.

EUROPEAN REPRESENTATIVE:



EMERGO EUROPE
Molenstraat 15
2513 BH, The Hague
The Netherlands
Phone: +31.70.345.6570
Fax: +31.70.346.7299

UNITED SURGICAL
2235 Pennsylvania St.
Fort Wayne, IN 46835
USA

United Ortho Enterprises

2235 Pennsylvania St. Fort Wayne, IN 46803
Phone: 800-227-8748; Fax: 877-688-7975
www.United-Ortho.com